

Announcing RMO's Exclusive

# Zero Corrector™

in collaboration with Dr. Peri Colino




## Introducing a Clinically Efficient and Versatile Class II and Class III Corrector Compatible with Both Aligner and Fixed Appliances

**Strategically Offset Bar Locates Forces Closer to Center of Resistance  
for Improved Bodily Movement and Reduced Tipping of Canines**

**Contoured Comfort  
Hook** for Easy Elastic  
Engagement

**Integrated Hook**  
Expands Treatment  
Options such as  
anterior space  
closure

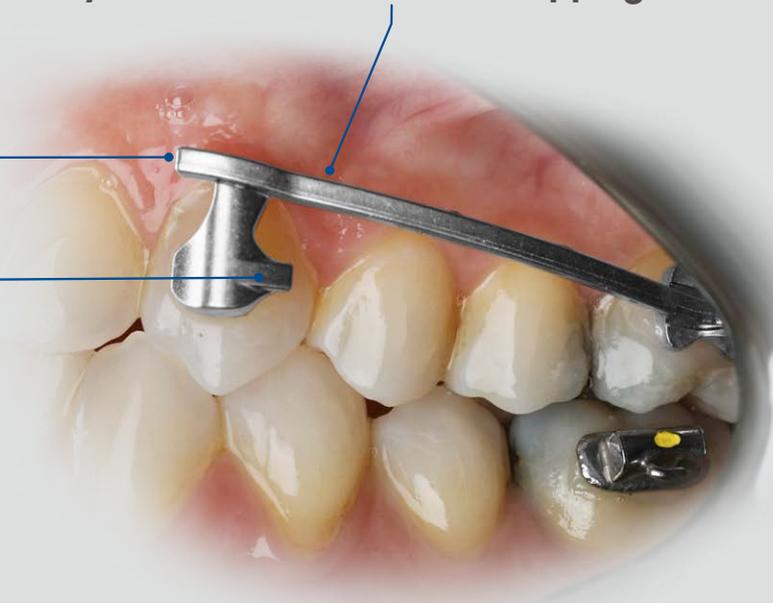

**Thoughtful Design Enhancements Ensure Optimal Performance,  
Ease of Placement and Patient Comfort**

**Rigid Yet Malleable  
Low Profile Arm**  
allows customizations

**Articulating Molar Pad** enhances  
molar rotating and uprighting while  
preventing over rotation

**Anatomical Cuspid Pad** ensures  
ideal fit for reliable bond strength

**Butterfly Shape Molar Pad** allows for  
precise positioning, easy removal of  
composite flash to avoid gingival irritation



# Class II Occlusal Correction



# Class III Occlusal Correction



## 6 Accurate Sizes Patient Sets – 2 Piece

| Size* | UL/LR, UR/LL |
|-------|--------------|
| 18mm  | Z02018       |
| 20mm  | Z02020       |
| 23mm  | Z02023       |
| 25mm  | Z02025       |
| 27mm  | Z02027       |
| 29mm  | Z02029       |



## Individual 1 Piece Per Package

| Size* | UL/LR  | UR/LL  |
|-------|--------|--------|
| 18mm  | Z01118 | Z01218 |
| 20mm  | Z01120 | Z01220 |
| 23mm  | Z01123 | Z01223 |
| 25mm  | Z01125 | Z01225 |
| 27mm  | Z01127 | Z01227 |
| 29mm  | Z01129 | Z01229 |



## Intro Kit – 6 Patient Sets

1 set of each size KZ1000  
 Contains 1 set each: 18mm, 20mm, 23mm, 25mm, 27mm, 29mm  
 2x Selection Guides

**\*Note:** Sizing is exact from long axis cuspid to first molar.

5x Scale Model of Zero™ Corrector for Patient Consultation

7 1/16" (179.4 mm) x 1 7/16" (36.5 mm)

REF Number

P00050



Typodont Zero™ Corrector

REF Number

P00031



Zero is a trademark of OrthoAmerica Holdings, LLC. © 2025 OrthoAmerica Holdings, LLC.



RMOOrtho.com  
 CustomerService@RMOOrtho.com  
 2165 Earlywood Dr. | Franklin, IN 46131

Toll Free: 800.525.6375  
 Fax: 303.592.8200  
 Local: 303.592.8209

Scan here to order your favorite RMO® products!  
 Order at [RMOOrtho.com](http://RMOOrtho.com)

